
Operation Manual

Laser Shutters

SR475 and SR476

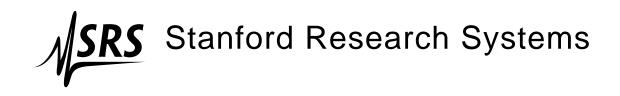

Revision 1.2 • September 1, 2011

Certification

Stanford Research Systems certifies that this product met its published specifications at the time of shipment.

Warranty

This Stanford Research Systems product is warranted against defects in materials and workmanship for a period of one (1) year from the date of shipment.

Service

For warranty service or repair, this product must be returned to a Stanford Research Systems authorized service facility. Contact Stanford Research Systems or an authorized representative before returning this product for repair.

Information in this document is subject to change without notice.

Copyright © Stanford Research Systems, Inc., 2011. All rights reserved.

Stanford Research Systems, Inc. 1290–D Reamwood Avenue Sunnyvale, CA 94089 USA Phone: (408) 744-9040 • Fax: (408) 744-9049 www.thinkSRS.com • e-mail: info@thinkSRS.com

Printed in U.S.A.

Document number 9-01636-903



Contents

G	enera	1 Information	iii
	Serv	rice	iii
	Env	ironment	iii
	Lase	er Safety	iii
	Sym	bols	iii
	Not	ation	iv
	Spee	cifications	\mathbf{V}
1	Ope	eration 1	-1
	1.1	Introduction to the Instrument	1-2
	1.2	Mounting	1-6
	1.3	Optical alignment	1 - 8
	1.4	Optical power capacity	1-9
	1.5	Operating environment	1-9
	1.6		1-10
	1.7	Shutter interface connector	1-12
2 Remote Operation		note Operation 2	-1
	2.1	Index of Commands	2-2
	2.2	Introduction	
	2.3	Commands	2 - 4
	2.4	Command List	

General Information

	The SR475 and SR476 Laser Shutters are high performance instru- ments designed to provide fast, versatile, and reliable mechanical blocking of a free-space laser beam while minimizing vibration and acoustic noise.
Service	
	Do not install substitute parts or perform any unauthorized modi- fications to this instrument. Contact the factory for instructions on how to return the instrument for authorized service and adjustment.
	Each SR475 and SR476 contains small parts and it is recommended that users do not remove the covers, and in particular, the sealed inner cover. There are no user-serviceable parts inside and reassembling the unit to its calibrated state may be difficult or impossible.
Environment	
	This product is intended for use only in a clean and dry laboratory environment. Operation in other environments may cause damage to the product.
Laser Safety	
	Certain hazards are always present when working with laser radia- tion. Visible and invisible beams of light have the potential to cause serious bodily injury including blindness or death and to cause sig- nificant damage to property.
MARNING	While the SR475 and SR476 are designed for use with laser systems, it is fully the users responsibility to ensure that safe operating conditions are maintained and to provide for fail-safe operation whenever an equipment failure could lead to a hazardous situation.

Symbols you may Find on SRS Products

Symbol	Description
\sim	Alternating current
	Caution - risk of electric shock
\rightarrow	Frame or chassis terminal
	Caution - refer to accompanying documents
↓-	Earth (ground) terminal
	Battery
\sim	Fuse
	On (supply)
0	Off (supply)

Notation

A

The following notation will be used throughout this manual.

A warning means that injury or death is possible if the instructions WARNING are not obeyed.

CAUTION

A caution means that damage to the instrument or other equipment is possible.

Typesetting conventions used in this manual are:

- Remote command names are set as A.
- Literal text other than command names is set as OFF.

Remote command examples will all be set in monospaced font. In these examples, data sent by the host computer to the SR475 and SR476 are set as straight teletype font, while responses received by the host computer from the SR475 and SR476 are set as *slanted* teletype font.



Specifications

Operation

General		
	Primary function	Block/unblock laser beam
	Controller/Driver type	SR470 Series
	Clear aperture	0.120" (3 mm) min. (SR475)
		0.040" (1 mm) min. (SR476)
	Lifetime	10^7 cycles min.
Shutter blade		
	Blade material	BeCu alloy, black oxide finish
	Optical power	Low-power laser beams; no absolute specification is established
	Bearing type	Sapphire jewel bearing
	Blade position	Servo controlled; neither open nor closed is preferred
	Power-off state	Blade position is indeterminate
	Hold time	Indefinite in either open or closed state.
Mounting		
0	Orientation	Any
	Mounting holes	8-32 (3), 1/4-20 (1)
	Housing material	Anodized aluminum
Performance		
Tenomance	Open/close time	5 ms full-stop to full-stop, in full-speed mode; User selectable to 10 ms, 20 ms, or 40 ms (SR475)
		4 ms full-stop to full-stop, in full-speed mode; User selectable to 8 ms, 16 ms, or 32 ms (SR476)
	Maximum rate	100 Hz (SR475) 125 Hz (SR476) (Limited by open/close time).
	Laser beam rise/fall time	Determined by beam size/quality and shutter mode. 500 μ s typ., in full-speed mode (SR475) 750 μ s typ., in full-speed mode (SR476)
	Insertion delay jitter	10 μs rms typ. (Measured at 10 Hz rep rate in full-speed mode.)
	Bounce	None allowed, at opening or closing

Electrical and Mechanical

General

Interface	SRS Shutter interface connector
Weight	0.28 lbs
Dimensions	$2.25'' \text{ W} \times 1.6'' \text{ H} \times 1.0'' \text{ D}$
Operating temperature	0 °C to 35 °C, non-condensing
Maximum temperature	90 °C

Shutter interface

Voltage inputs	+12 VDC, +4.5 VDC
Current draw	1.20 A max. continuous (12 VDC)
	0.25 A max. continuous (4.5 VDC)
Signal lines	Serial TX, Serial RX, logic-level control
Connector type	Mechanically compatible
	with 6-pin IEEE-1394 cables.
Cable length	10 ft (3 m) max.
Cable shield	Internally connected to housing
	•





1 Getting Started

This chapter gives you the necessary information to begin using the SR475 and SR476 Laser Shutters.

In This Chapter

1.1	Introd	duction to the Instrument
	1.1.1	Overview 1-2
	1.1.2	Design and theory of operation $\ldots \ldots \ldots \ldots 1-2$
	1.1.3	Mechanical performance
	1.1.4	Power-off and standby states
	1.1.5	Power-on state 1–6
	1.1.6	Resetting the shutter head
1.2	Moun	nting 1-6
	1.2.1	Vibration considerations
	1.2.2	Thermal considerations 1–7
	1.2.3	Grounding considerations
1.3	Optic	al alignment 1-8
1.4		al power capacity 1–9
1.5		ating environment
	1.5.1	Dust 1-9
	1.5.2	Magnetic field
	1.5.3	Temperature
	1.5.4	Vacuum
	1.5.5	Oxygen-rich environments
1.6	Alterr	nate speed modes
1.7	Shutt	er interface connector
	1.7.1	Direct interfacing

1.1 Introduction to the Instrument

The SR475 and SR476 Laser Shutters are high performance instruments designed to provide fast, versatile, and reliable mechanical blocking of a free-space laser beam while minimizing vibration and acoustic noise.

The shutters are suitable for low-power beam switching and pulse selection applications. They are also particularly suited to provide additional extinction to optical systems incorporating fast optoelectronic shutters.

1.1.1 Overview

The basic function of the laser shutter, or "shutter head" is to block or unblock a laser beam when commanded to do so.

The shutter head is normally operated in conjunction with a shutter driver or controller that provides power, a serial communications interface, and a control input. Details about this interface are given in §1.7.

1.1.2 Design and theory of operation

Mechanical laser shutters have traditionally been electronic descendants of camera shutters, consisting of a solenoid-driven blade (or set of blades) with a return spring. The shutter's open and closed positions are then determined by physical stops that limit the blade travel.

By contrast, the SR475 and SR476 Laser Shutters do not use physical stops to control the blade position. Instead, the blade position is precisely measured in real time and a digital servo loop is used to hold the blade in the desired location. During a transition between the open and closed states, the shutter blade is guided along a predetermined trajectory such that it arrives at its final position with zero velocity, without encountering a physical stop.

This design reduces or eliminates many of the less desirable characteristics of conventional solenoid-based mechanical laser shutters. Most importantly, it removes the high impulse that is encountered when a solenoid or a return spring drives a shutter blade into a hard stop. This high impulse is responsible for the bulk of vibration and mechanical noise associated with shutters, but is also a contributing factor in the limited lifetime of such devices.

A typical shutter with a solenoid and return spring necessarily has a defined *normal* state—either open or closed—and expends energy not only to reach but also to hold in the normal state. By using an

active servo system instead of a solenoid and return spring, the SR475 and SR476 do not have a preferred state, and can stay indefinitely in either state without significant power dissipation. Significant power is only used when opening or closing, and each shutter head is rated for continuous duty without duty cycle restrictions.

The shutter blade in an SR475 or SR476 is mounted in a sapphire jewel bearing to minimize friction, and is driven as a part of a tightly integrated rotary motor. The motor design, along with the impulse reduction from eliminating physical stops, results in a head lifetime that is orders of magnitude more than is typically found in other laser shutter systems.

The on-board electronics in the shutter head detect the blade position, servo control the blade position, and provide control, communication, and error detection functions. Electrical, over-temperature, and mechanical faults are detected and can be reported over the shutter head's serial interface.

1.1.3 Mechanical performance

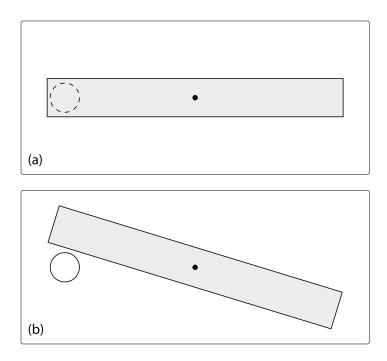


Figure 1.1: Mechanical action of the SR475 and SR476 shutter heads. Part (a) illustrates the shutter head in the *closed* position, where the shutter blade blocks the aperture. Part (b) shows the shutter head in the *open* position, where the shutter blade is clear of the aperture.

The location of the shutter blade in its open and closed positions is illustrated in Fig. 1.1. In the open position the shutter blade is clear of

the laser aperture, and in the closed position the shutter blade fully covers the aperture.

The process of opening or closing the shutter consists of guiding the blade along a smooth, pre-determined trajectory between the closed and open positions. In what follows we will describe this in more detail. For clarity we will discuss only the case of opening, but the sequence of events applies equally to closing the shutter.

When the shutter is initially closed and it receives a signal to open, it begins moving the shutter blade under servo-loop control along the preprogrammed "opening" trajectory. Once this has begun, the shutter will guide the blade along the entire trajectory—bringing the blade to a full stop at the open position—before carrying out any subsequent instruction to close. In other words, a transition in progress cannot be reversed or aborted (except by disabling the motor).

Sample trajectories for opening and closing are illustrated in Fig. 1.2. For the SR475, each pre-programmed trajectory takes 5 ms to complete when running at full speed. This means that the minimum time to complete a closed-open-closed full cycle is 10 ms, which gives the shutter's maximum rate of 100 Hz. The minimum exposure time between the midpoint of opening and the midpoint of closing is 5 ms. For the SR476, each full-speed pre-programmed trajectory takes 4 ms, and a full cycle takes a minimum of 8 ms. This gives a maximum rate of 125 Hz and a minimum exposure time of 4 ms. For both models, when full speed is not needed, there are alternate speed modes available that use less power and create less vibration; see §1.6.

1.1.4 Power-off and standby states

The SR475 and SR476 Laser Shutters do not incorporate return springs. The shutter blade is in an indeterminate position whenever the shutter is powered off or in the *Standby* mode.

The shutter is in the *Standby* mode whenever its +4.5 VDC "logic" power supply is available but the motor is not enabled. In this mode, the logic power supply provides power for serial communications, but the motor is turned off and the shutter blade position is therefore indeterminate.

The shutter head can be placed into *Standby* with a remote command. It also enters *Standby* mode automatically when certain fault conditions are detected. These conditions include mechanical failure, over temperature, and undervoltage on the +12 VDC power supply for the motor. The shutter head may also be put into *Standby* by intentionally turning off the +12 VDC power supply. To ensure that

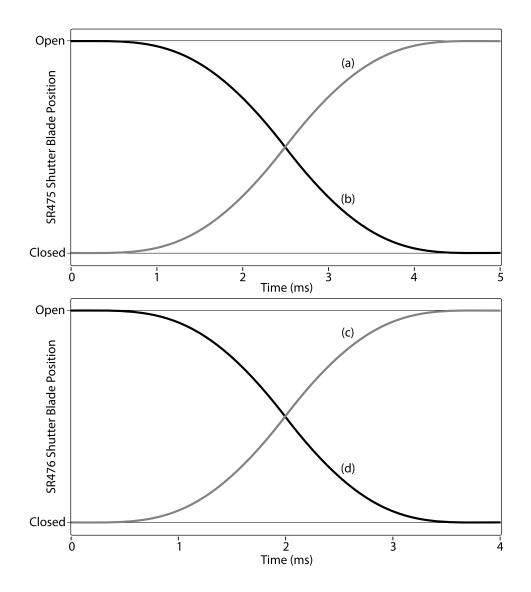


Figure 1.2: Blade trajectory for the SR475 Laser Shutter (top) and SR476 Laser Shutter (bottom). Curve (a) illustrates the position of the SR475 blade as it moves from the closed to the open position, and curve (b) illustrates the closing trajectory. Curve (c) illustrates the position of the SR476 blade as it moves from the closed to the open position, and curve (d) illustrates the closing trajectory.

the shutter head enters standby mode, the +12 VDC power should be kept low for at least three seconds. Once in *Standby* mode, the shutter head will remain so until reset.

Note: The +4.5 VDC "logic" power supply is always provided to a shutter while it is connected to a powered-on SRS shutter controller such as the SR470 Shutter Controller or the SR474 Four-Channel Shutter Driver. Placing a shutter drive channel on one of these instruments into the "off" or "sleep" states will put a SR475 or SR476 attached to that channel into the *Standby* mode.

1.1.5 Power-on state

At turn-on, the SR475 or SR476 executes a brief mechanical initialization to verify motor function and perform the initial servo-lock of the shutter blade position. This can potentially last several seconds and/or involve minor mechanical bumps in excess of normal operation.

All settings and modes of the SR475 and SR476 are in volatile memory, and at power-on the instrument always returns to its default state.

1.1.6 Resetting the shutter head

The SR475 and SR476 can be reset by turning off and restoring the +4.5 VDC power supply, through a remote command to reset the shutter head, or by restoring the +12 VDC power supply after it has been turned off.

1.2 Mounting

SR475 and SR476 Laser Shutters can be mounted and operated in any orientation. 8-32 mounting holes are provided on three sides of the shutter head, and a single 1/4 - 20 mounting hole is provided on the bottom side, as shown in Fig. 1.3.

Orientation may be chosen for convenient routing of the shutter cable, for aligning to the appropriate beam height above an optical table, or for pulsing one of two closely separated laser beams.

1.2.1 Vibration considerations

SR475 and SR476 Laser Shutters should be firmly mounted to a solid surface. Using long and/or narrow mounting posts will increase the apparent vibration induced by the shutter.

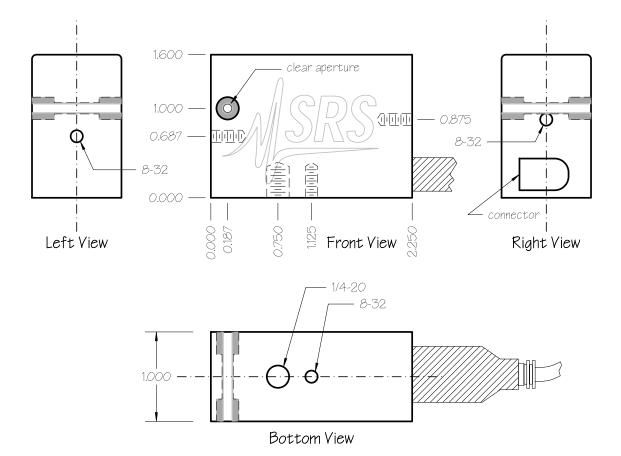


Figure 1.3: Physical dimensions and mounting of the SR475 and SR476. The clear aperture is 0.120" (3 mm) for the SR475 and 0.040" (1 mm) for the SR476.

For optimal results, use short, large-diameter mounting posts or other sturdy optomechanical hardware to mount SR475 and SR476 Laser Shutters to appropriate surfaces. Attaching the shutter head rigidly to a heavy, damped optical table will minimize its vibrational effects.

1.2.2 Thermal considerations

Like most other mechanical shutters, the performance of the SR475 and SR476 Laser Shutters are fundamentally limited by their ability to dissipate heat. The motor can get hot when run continuously (or nearly so), particularly in full-speed mode.

The outer case of the shutter head may get warm but should not become too hot to touch under normal operating conditions. Solidly mounting shutter heads to metal surfaces can assist heat sinking



and help them to run at a comfortable temperature, particularly in a warm laboratory environment.

An internal thermal sensor will trip off if the shutter head overheats. In this case, the head will declare a fault and go to the *Standby* state.

1.2.3 Grounding considerations

For reliable operation of a SR475 or SR476, it is essential that the shield of the shutter cable is properly grounded.

Inside the shutter head, there is an electrical connection between its chassis (case) and the shield of the cable. If the chassis of the shutter head is electrically connected to earth ground, for example by way of attachment to a grounded optical table, then the shield will be grounded as well.

SRS shutter controllers, such as the SR470 Shutter Controller and the SR474 Four-Channel Shutter Driver, provide optional internal grounding of the shutter cable shield. If a shutter case is properly grounded and thereby providing a suitable ground to the cable shield, it may be desirable to disconnect the shield-ground connection within the controller; consult that instrument's operation manual for details.

Apart from the cable shield, the shutter cable uses a separate "signal ground" as one of the power supply lines. This ground is normally isolated from the shutter case and cable shield. Under certain circumstances it may be desirable (for the control of ground loops, for example) to establish an electrical connection between the signal ground and the cable shield. To do so, unplug the shutter head and solder a jumper wire across location J11 on the circuit board. To access the circuit board, you will need to temporarily remove the outer case, which is held on by four screws.

1.3 Optical alignment

The SR475 and SR476 Laser Shutters are designed for easy alignment. The clear aperture of the SR475 is aperture of 0.120" (3 mm), and the clear aperture of the SR476 is 0.040" (1 mm). The laser beam may enter the shutter head from either side, but there is less potential for leakage of indirect scattered light if the beam enters from the side marked "Front" in Fig. 1.3.

SRS shutter controllers incorporate an "Align" mode that chops the shutter open and closed at a rate of 1 Hz to assist with manual alignment. A similar chopping function is also directly built into the shutter head firmware, and can be enabled over the remote interface.

The rise and fall time of a laser beam chopped by the shutter are application dependent, determined by the size and profile of the laser beam. A narrower beam will blink on and off quickly compared to a broader beam, since the blade edge will take less time to transit a smaller beam. Similarly, the insertion delay, i.e., the time between initiating an open-closed transition and when the laser beam is actually switched, depends upon the beam profile and its location in the clear aperture. In general, the insertion delay will be different for opening and closing.

While it may be desirable to gently focus a laser beam through the shutter aperture to reduce the apparent open/close time, it should be noted that this speed increase comes as a tradeoff with increased sensitivity to insertion delay jitter.

1.4 Optical power capacity

The SR475 and SR476 use a thin BeCu alloy shutter blade with a black oxide finish and are intended for use in "low-power" laser beam applications.

The maximum power handling capability of a laser shutter depends on wavelength, duty cycle, and many other environmental factors. SRS does not establish optical power ratings for laser shutters; please contact SRS if you would like to request sample shutter blades to evaluate for your application.

1.5 Operating environment

The SR475 and SR476 Laser Shutters are designed for operation in a clean and dry laboratory environment.

1.5.1 Dust

Both the bearing and the position detection system within the shutter head are susceptible to impairment in a dusty environment. This can potentially result in a fault condition, when either the blade is unable to move to the commanded position, or when the readout system is unable to verify that the blade is in the correct position.

Always operate the laser shutter with its outer cover in place to minimize the potential for dust to collect inside the shutter head.

If dust contamination inside the shutter is suspected, the outer cover can be unscrewed and dust can be *gently* blown out. Use only lowvelocity air to remove dust from the shutter head, for example a handheld hair dryer with the heat turned off.



npressed otentially
echanism ; external
ation in a
exceeds ne circuit <i>Standby</i>
r head is lt.
use in a onvection l cooling

Also note that adhesives used in construction of the SR475 and SR476 Laser Shutters are not rated for vacuum compatibility.

1.5.5 Oxygen-rich environments

Because of the adhesives used in construction, the SR475 and SR476 Laser Shutters are not suitable for use in oxygen-enriched, pure oxygen, or other similar environments.

1.6 Alternate speed modes

The SR475 and SR476 Laser Shutters are by default configured to operate in *full-speed mode*, where the opening/closing time is fixed at its minimum possible value. During the transition from one state to the other, the shutter blade is guided along a pre-determined trajectory between the two states, as discussed in §1.1.3.



Three alternate modes of operation are provided, each of which guides the shutter along the same trajectory but at a slower rate. This leads to a significant reduction in power usage and physical impulse. These slower modes can be used in many applications where absolute speed is not as crucial as minimizing vibration.

Mode	Transition time	Max Rate
0	5 ms	100 Hz
1	10 ms	50 Hz
2	20 ms	25 Hz
3	40 ms	12.5 Hz

Table 1.1: Model SR475 Laser Shutter speed modes

Mode	Transition time	Max Rate
0	4 ms	125 Hz
1	8 ms	62.5 Hz
2	16 ms	31.25 Hz
3	32 ms	15.625 Hz

Table 1.2: Model SR476 Laser Shutter speed modes

The actual transition times in the different modes for the SR475 and SR476 Laser Shutters are given in Table 1.1 and Table 1.2, respectively. Mode **0** is the full-speed mode, and modes 1–3 are alternate modes. The corresponding maximum rates are for a full cycle including an opening and closing transition.

The alternate modes can be selected using the remote interface; please see §2.4 for more information.



1.7 Shutter interface connector

The SR475 and SR476 Laser Shutters are normally operated through SR470 Series Laser Shutter drivers (such as the SR470 Shutter Controller or the SR474 Four-Channel Shutter Driver) that provide power and external interfaces. Users may also wish to directly interface the SR475 or SR476 to their own systems without the use of additional hardware.

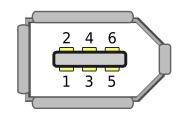


Figure 1.4: Shutter connector pinout, as seen looking into the connector.

Connector Pin	Description
1	+12 VDC in (motor power supply)
2	Signal ground / power return
3	+4.5 VDC in (logic power supply)
4	Control signal in (logic level)
5	RX (Serial input)
6	TX (Serial output)

Table 1.3: Shutter connector pinouts

The pinout for the connector is shown in Fig. 1.4 and the signal names are shown in Table 1.3.

The shield of the shutter connector is connected to the case of the shutter head. Please see §1.2.3 about grounding considerations for more information.

1.7.1 Direct interfacing

 The SR475 and SR476 have limited or no internal protection against reverse polarity, missing power, or overvoltage on the power supply pins. The misapplication of power may cause circuit damage. For most applications, SRS recommends using the SR475 or SR476 together with an SR470 Series Laser Shutter driver such as the SR470 Shutter Controller or the SR474 Four-Channel Shutter Driver.

The shutter connector is mechanically compatible with 6-pin IEEE-1394 cables. Clean, well-regulated supply voltages of +4.5,+12 VDC must be provided at the shutter head following the pin-out specified in Table 1.3. Ground must be provided on pin 2, with chassis ground on the shield.

Note that standard IEEE-1394 cables are crossover cables that reverse the signals of pins (3,5) and of (4,6) at the two ends.

1.7.1.1 Serial communication

The features and operation of the SR475 and SR476 remote interface are discussed in §2.

The shutter head serial interface accepts logic-level input (3.3 VDC, 5 VDC tolerant) signals on the RX pin and transmits signals at 3.3 VDC on the TX pin. The protocol is fixed at 19200 baud, with 8 data bits, 1 stop bit, no parity, and no flow control. The polarity of both signals is such that the start bit and logical "0" are at ground potential, while the logical "1" is at= +3.3 VDC.

If the user intends to directly establish serial communication with the shutter head, independent of an SRS shutter controller, it may be necessary to use an external signal level converter to translate between the logic-level signals at the shutter head and (for example) the voltages and "inverse" polarity of an RS-232 interface.

1.7.1.2 Control signal

The logic-level control signal is the most important input to the shutter head. Unless overridden through the remote interface, a logical high signal on this input commands the shutter to the open state, while a low signal commands the shutter to the closed state.

When a change in the input signal is detected, the shutter will immediately begin a transition to the new state if a transition is not already in progress. A pending state change request will be executed only after completion of the transition. Changing the level of the control signal introduces the lowest possible jitter (phase noise) of any method of initiating a transition, and is the method used by SR470 Series Laser Shutter drivers.

If directly interfacing to the SR475 or SR476, the control signal accepts a logic-level input (3.3 VDC, 5 VDC tolerant). To ensure reliable operation, control pulses to the SR475 or SR476 should be at least 200 μ s in duration.



2 Remote Operation

This chapter describes operating the SR475 and SR476 over the builtin serial interface.

In This Chapter

2.1	Index	of Commands	2 - 2
2.2	Introc	luction	2-3
	2.2.1	Control modes	2-3
	2.2.2	Timing considerations	2-3
2.3	Comn	nands	2-4
	2.3.1	Command syntax	2 - 4
2.4	Command List		
	2.4.1	Shutter motion control commands	2 - 5
	2.4.2	Shutter control mode commands	2 - 6
	2.4.3	Communication control commands	2 - 7
	2.4.4	Utility commands	2 - 8
	2.4.5	Shutter queries	2 - 8

2.1 Index of Commands

Motion control		
@	2 - 5	Open
A		Close
В	2-5	Toggle state
J	2 - 5	Toggle alignment mode
0	2-6	Toggle alignment mode Select speed mode 0
1	2-6	Select speed mode 1
2	2-6	Select speed mode 2
3	2-6	Select speed mode 3
Control mode	command	S
E	2-6	Serial control lockout
F	2-6	Logic-level control lockout
G		Flexible (default) control mode
Communicatio	on control	commands
L	2 – 7	Enable assertive mode
Μ	2 - 7	Disable assertive mode
Utility commai	nds	
Utility comma C		Reset
-	2 - 8	Reset Standby
C	2 - 8 2 - 8	
C K	2-8 $2-8$ $2-8$	Standby
С К О	2-8 2-8 2-8	Standby
C K O Shutter querie	2-8 2-8 2-8	Standby Assert fault State
C K O Shutter querie S	2 - 8 2 - 8 2 - 8 s 2 - 8 2 - 8 2 - 9 s 2 - 8 2 - 9 s 2 - 9 s 2 - 9 s 2 - 9 s 2 - 9 - 9 2 - 9 2 - 9 2 - 9 2 - 9 2 - 9 2 - 9 - 9	Standby Assert fault State Rate
C K O Shutter querie S R	2 - 8 2 - 8 2 - 8 2 - 8 s 2 - 8 2 - 9	Standby Assert fault State
C K O Shutter querie S R T	2 - 82 - 82 - 83s $2 - 82 - 92 - 92 - 92 - 9$	Standby Assert fault State Rate Temperature
C K O Shutter querie S R T W	2 - 82 - 82 - 8s $2 - 82 - 92 - 92 - 92 - 92 - 92 - 10$	Standby Assert fault State Rate Temperature Error Word

2.2 Introduction

When an SR475 or SR476 is operated through an SRS shutter controller such as the SR470 Shutter Controller or the SR474 Four-Channel Shutter Driver, that instrument normally signals the shutter head to open or close through the shutter head's logic-level control input.

The SR475 and SR476 also support serial communication and may be remotely programmed via this interface.

Any host computer interfaced to an SR475 or SR476 shutter head through a shutter controller may use that controller to relay (pass through) serial commands to the shutter head. This allows the SR475 or SR476 to be controlled remotely through the controller's interfaces, which may include Ethernet, RS-232, and/or GPIB.

The SR475 and SR476 can also be operated as stand-alone devices without an SRS controller. See §1.7.1.1 for details.

2.2.1 Control modes

By default, the SR475 and SR476 follow the state dictated by the logic-level control input (§1.7.1.2). However, this behavior can be overriden through the use of the serial interface.

Executing any remote command that changes the position of the shutter blade will force the shutter into *logic level lockout mode*, where it will no longer respond to the logic-level control signal. Note that SRS shutter controllers, such as the SR470 Shutter Controller and the SR474 Four-Channel Shutter Driver, normally control the state of the shutter through the logic-level input, but cannot do so when that interface is locked out. Hence, care should be exercised when directly commanding the shutter state over the serial interface.

Logic level and remote interface control of the blade position can be separately enabled or disabled through the remote interface.

Also keep in mind that settings and modes of the SR475 and SR476 only exist in volatile memory, and at power-on the instrument always returns to its default state.

2.2.2 Timing considerations

While it is possible to open and close the SR475 and SR476 Laser Shutters through direct commands over its serial interface, these operations are executed with significantly lower timing precision than control through the logic-level input. Whenever timing is important, serial control of the shutter position should be avoided.

<u>SRS</u>

2.3 Commands

This section provides syntax and operational descriptions for remote commands.

2.3.1 Command syntax

Communication with the SR475 and SR476 uses ASCII characters. All commands are one character long and are case sensitive.

Commands are executed when the command character is received without waiting for additional data or confirmation. No command terminator should be used. While this style of execution is referred to as "immediate," remote commands are not executed with precision timing (§2.2.2).

The command buffer is limited to 16 bytes. If the buffer overflows, a bit will be set in the error word, and the last command will be discarded.

Query commands, which are commands that request a return of data, will return exactly seven bytes of data, of which the last byte is $\langle LF \rangle$ (ASCII 10). Numbers are formatted as signed integers, where the leading character is '-' for negative numbers or ' ' (ASCII 32) for positive numbers.

Aside from communication errors, commands may fail due to either syntax or execution errors. Not all such failures will generate a fault. As an example, if the shutter is initially open, a command to open the shutter will simply be ignored.

All commands (ASCII characters) that are not in the command list are reserved for future use and/or factory use. Commands that are not in the list may produce unpredictable results and should be avoided.

Errors and changes in shutter status can be detected by requesting the status word and error word, with commands W and Z.

2.3.1.1 Examples

Commands may be given with a simple example illustrating their usage. In these examples, commands sent by the host computer to the SR475 or SR476 are set as straight teletype font, while responses received the host computer from the shutter head are set as *slanted teletype font*.

These examples are not exhaustive, but are intended to provide a convenient starting point for user programming.

2.4 Command List

2.4.1 Shutter motion control commands

@	Open
	Moves the shutter blade from the closed state to open state.
	When executed, this places the shutter into logic level lockout mode. (See description of command F.)
	The command to open the shutter will be ignored unless (1) the shutter is in the closed position, and (2) serial control of the blade position has not been disabled with command E.
A	Close
	Moves the shutter blade from the open state to closed state.
	When executed, this places the shutter into logic level lockout mode. (See description of command F.)
	The command to close the shutter will be ignored unless (1) the shutter is in the open position, and (2) serial control of the blade position has not been disabled with command E.
В	Toggle state
	Moves the shutter from the Open state to Closed state or vice versa, depending on the initial state.
	When executed, this places the shutter into logic level lockout mode. (See description of command F.)
	The command to toggle the shutter state will be ignored unless (1) the shutter is in the open or closed position, and (2) serial control of the blade position has not been disabled with command E.
J	Toggle alignment mode
	Enable alignment mode if it is presently disabled. Disable alignment mode if it is presently enabled.
	In alignment mode, the shutter continuously opens and closes at fixed rate of approximately 2 Hz. This mode can be helpful for optical alignment of the shutter head; see §1.3.
	When executed, this places the shutter into logic level lockout mode. (See description of command F.)
	This command will be ignored unless (1) The shutter is either in the open or closed position or already in alignment mode, and (2) serial



	control of the blade position has not been disabled with command E.
0	Select speed mode 0
	Select full-speed mode for the shutter head. This is the default setting.
	The transition speed (and corresponding mechanical vibration) can be changed using mode commands 0, 1, 2, and 3. The actual speeds for these modes, for the SR475 and SR476 Laser Shutters, are given in §1.6.
1	Select speed mode 1
	Select alternate speed mode 1 for the shutter head, so that it moves at one-half of full speed between open and closed.
2	Select speed mode 2
	Select alternate speed mode 2 for the shutter head, so that it moves at one-quarter of full speed between open and closed.
3	Select speed mode 3
	Select alternate speed mode 3 for the shutter head, so that it moves at one-eighth of full speed between open and closed.

2.4.2 Shutter control mode commands

E	Serial control lockout
	Disable remote (serial command) control of shutter blade position. Commands that do not change the blade position are not affected.
	When this mode is selected, the blade position can only be controlled through the shutter head's logic-level control input. Serial control can be re-established through commands F, G, or through a hard reset with command C or by removing power (§1.1.5).
F	Logic-level control lockout
	Disable logic-level control of shutter blade position.
	When this mode is selected, the blade position can only be controlled through the remote (serial) interface. Logic-level control can be re- established through commands E, G, or a hard reset with command C or by removing power (§1.1.5).

G	Flexible (default) control mode		
	Serial control lockout and logic-level control lockout are both dis- abled, and either may control the shutter blade position.		
	This is the default control mode, and is always selected after the shutter head is reset.		
2.4.3 Communic	ion control commands		
L	Enable assertive mode		
	Places the shutter in assertive communication mode.		
	Under normal operation, the shutter head will never initiate com- munication over its serial interface; it will only respond to queries.		
	When assertive mode is enabled, the shutter head will actively re- quest communication—by issuing a Serial Break character—if a "fa- tal" fault occurs. The fatal faults are those that will trip the shutter into the standby state: temperature, 12 V, position, and motor fail- ures. A fatal fault can also be simulated by issuing command O.		
	While this mode is enabled, all queries that request data are ignored. This mode is disabled by either issuing command M or by a hard reset, either with command C or by removing power (§1.1.5).		
M	Disable assertive mode		
	Disable assertive communication mode.		
	If assertive mode is currently enabled, this command will turn it off and restore the ability to communicate normally with the shutter head.		



2.4.4 Utility commands

С	Reset		
	Performs a hard reboot of the shutter head. Upon restart the head will initialize and attempt to move the blade to the state indicated by the logic-level control input. The full startup sequence may take up to one second.		
K	Standby		
	Turns off the shutter motor and places the shutter head in the standby state, as described in §1.1.4. The shutter blade position becomes indeterminate. Once this mode is entered, the blade position will remain indeterminate and unpowered until a hard reset, either with command C or by removing power (§1.1.5).		
0	Assert fault		
	This command simulates a "fatal" fault which will trip the shutter to the standby mode. This function is provided as a utility for reliability testing.		
2.4.5 Shutter queries			
S	State		
	Query the current position of the shutter blade. The returned value may indicate that the shutter is open, closed, or indeterminate. The shutter is said to be in an indeterminate position (i.e., not open or closed) while in transit between open and closed, or when powered off in standby mode (§1.1.4).		
	ResultMeaning1Open0Closed-1Indeterminate		
Example:	S 1 The returned value of 1 indicates that the shutter is presently open.		

R		Rate		
		Query the maximum rate of the shutter head.		
		The returned value indicates the maximum rep rate of the shutter head, in Hz, rounded down to the nearest integer. The rate depends on the speed mode that is selected; see §1.6.		
	Example:	R 100 The shutter has a maximum rate of 100 Hz.		
	Example:	R 12 This value is returned when the shutter is in alternate speed mode (mode 3), which has a maximum rate of 12.5 Hz.		
Т		Temperature		
		Query the shutter temperature.		
		The returned value is the approximate temperature as measured at the printed circuit board in the shutter head, in °C.		
	Example:	T		
	-	42 The shutter is at an approximate temperature of 42 °C.		
W		Error Word		
		Query the shutter head error word.		
		The error word is a 15-bit value that is returned as an unsigned integer. This query also clears the four communications error bits.		
		Bit Meaning		
		• Serial buffer overflow		
		1 Serial syntax error		
		2 Serial communications error		
		3 Serial Timeout4 (Reserved)		
		5 Temperature fault		
		6 12 V fault		
		7 Motor fault		
		8 Firmware fault		
		9 Position fault		
		10-15 (Reserved)		
		Bits 5, 6, 7, and 9 correspond to "fatal" faults that will trip the shutter to the standby state. Of the fatal faults, only the first one that is		



	Example:	W 64	will be recorded in the error word.	
		The retu power s	rned value of 64 indicates a problem with the 12 VDC shutter upply.	
Х		Model Number		
		5	he shutter head model. Returns the 7 character ASCII idenher " $SR475(LF)$ " or " $SR476(LF)$ ".	
	Example:	X SR475		
Y		Serial Number		
		Query th	ne shutter head serial number.	
	Example:	Y <i>1234</i> The retu 1234.	rned value indicates that the shutter head serial number is	
Z		Status W	Vord	
		an unsig	ne shutter head status word, a 15-bit value that is returned as med integer. It can be used to verify the current modes and of the shutter head.	
		Bit	Meaning	
		0	Shutter OK	
		1	Motor Enabled	
		2 3	Direction In transit	
		4	Logic-level lockout	
		5	Serial lockout	
		6	Alignment mode	
		7,8	(Reserved)	
		9	Assertive mode	
		10	(Reserved)	
		11 12	Locked Mode0 (Lo bit of mode number)	
		12	Model (Lo bit of mode number) Model (Hi bit of mode number)	
		14,15	(Reserved)	
			<i>utter OK</i> indicates that the shutter head is working properly.	
		DI. 4 3 5		

Bit 1, *Motor Enabled* is 1 when the shutter motor is enabled.

Bit 2, *Direction* is 1 when the shutter is currently commanded to the open state, and 0 when commanded to the closed state. (This bit is only meaningful when the motor is enabled.)

Bit 3, *In transit* is 1 when the shutter blade is between the open and closed positions. (This bit is only meaningful when the motor is enabled.)

Bit 4, *Logic-level lockout* is 1 when the logic-level control input is disabled.

Bit 5, *Serial lockout* is 1 when the serial control of the blade position is disabled.

Bit 6, *Alignment mode* is 1 when the shutter head is in alignment mode. (This bit is only meaningful when the motor is enabled.)

Bit 11, *Locked* is 1 when the shutter blade is servo locked to the correct position.

Bits 12 and 13, *Mode number* give a number from 0 to 3 that indicates which of the alternate speed modes is presently selected.

