
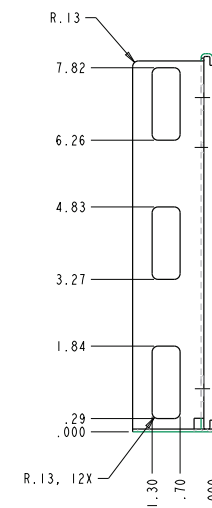


REVISIONS			
REV	DESCRIPTION	DATE	APPROVED
B	ADDED NOTCH	5/6/03	GAW
C	CHANGED MATERIAL AND FINISH	8/18/07	GAW
D	ADDED SILKSCREEN & CHANGED HOLE SIZES	2/27/08	GAW
D1	ADDED RoHS NOTE	10/6/17	GAW

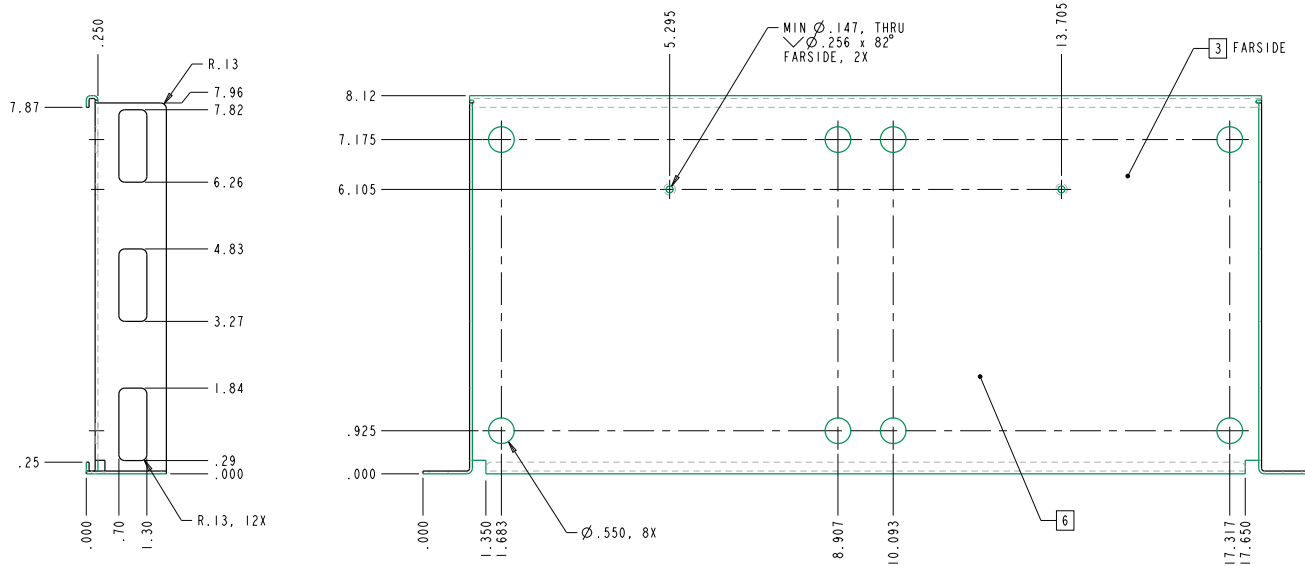

NOTES: UNLESS OTHERWISE SPECIFIED

- MATERIAL: 16 GA. CRS, (.059 THICK).
- FINISH: ZINC PLATE CLEAR PRIOR TO PAINT PAINT ENTIRE PART INCLUDING EDGES, WARM GREY CARDINAL #8101-50714, MEDIUM TEXTURE OVERSPRAY PERMISSIBLE.
- PART MARKING: VENDOR ID, MONTH/YEAR OF MANUFACTURE, PART NUMBER, AND CURRENT REVISION LEVEL USING METAL AND/OR RUBBER STAMP, APPROXIMATELY WHERE SHOWN.
- REMOVE ALL BURRS AND SHARP EDGES.
- BEND RADIUS SHALL BE .03 MAX.
- SILKSCREEN TOP USING SRS ARTWORK FILENAME: 7-01480 TOP SILKSCREEN.ai, BLACK EPOXY.
- PART TO BE RoHS COMPLIANT.

NEXT ASSY		SR445A	DO NOT SCALE DRAWING	23-Apr-03	SCALE: 0.75	SHEET: 1 OF 1
-----------	--	--------	----------------------	-----------	-------------	---------------

MATERIAL		UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES.		SRS STANFORD RESEARCH SYSTEMS, INC. 1290-D Redwood Avenue, Sunnyvale CA 94089	
FINISH		TOLERANCES: .X ± .02 .XX ± .010 ANGLE ±		DRAWING TITLE: TRAY, RACK MOUNT, SR445A	
ENG BY/DATE G. WATERS		SIZE/PROJ FILENAME: D		SRS PART NO. 7-01480-720	
OWN BY/DATE G. WATERS		23-Apr-03		REV D1	



CAUTION: SECURE UNIT
WITH 6-32 SCREW FROM BELOW



STANFORD RESEARCH SYSTEMS
FS730 DUAL RACK MOUNT

CAUTION: SECURE UNIT
WITH 6-32 SCREW FROM BELOW